How to use the Flare Water Flosser

1. Attach the Nozzle: Insert the nozzle into the top of the flosser until you hear a click. Replace the nozzle every 3 months or if it becomes loose.


2. Fill the Water Tank: Open the lid of the water tank and pour in water or your preferred solution.

3. Close the Water Tank: Securely close the water tank lid.

4. Position the Nozzle: Place the nozzle tip above the gum line between your teeth. Gently close your lips around the nozzle to prevent splashing.
5. Turn On the Flosser: Press the power button. The LED light will indicate that the flosser is ready. Select the desired mode (Soft, Standard or Pulse).

6. Floss Your Teeth: Move the nozzle along the gum line and between teeth. Ensure you clean all interdental spaces, including behind the back teeth.

7. Refill if Necessary: If needed, refill the water tank and continue flossing.
- Visit be-lucent.com for replaceable nozzle tips.

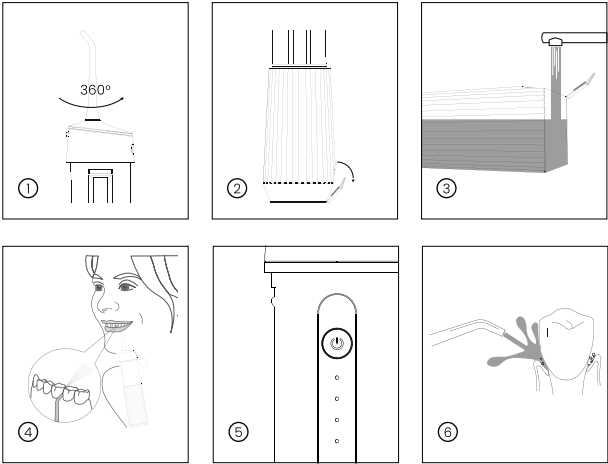LUCENT

FLARE

FLARE

Water Flosser

Instructions Of Use	English	5-6
Handbuch	Deutsch	7-8
Mode D'emploi	Français	9-10
Manual De Usuario	Español	11-12
Manuale D'uso	Italiano	13-14
Gebruiksaanwijzing	Nederlands	15-16
Bruksanvisning	Svenska	17-18
Brugsanvisning	Norsk	19-20
Brugsanvisning	Dansk	21-22
Käyttöohjeet	Suomalainen	23-24
Instrukcja Użycia	Polski	25-26

INTRODUCTION


Thank you for buying the Flare Water Flosser.

Flare is designed to provide the most convenient and efficient care for your oral hygiene. The push-button switch allows you to control the power and mode with ease. Just floss, smile, and keep shining.

Of course, the water flosser is waterproof, so you don't have to worry about getting it wet. Maybe even floss in the shower.

It's important to always use a fresh nozzle tip. Dentists recommend replacing it at least every 3 months. For best performance, we recommend changing as soon as the tip shows signs of wear.

Visit be-lucent.com for replaceable nozzle tips.

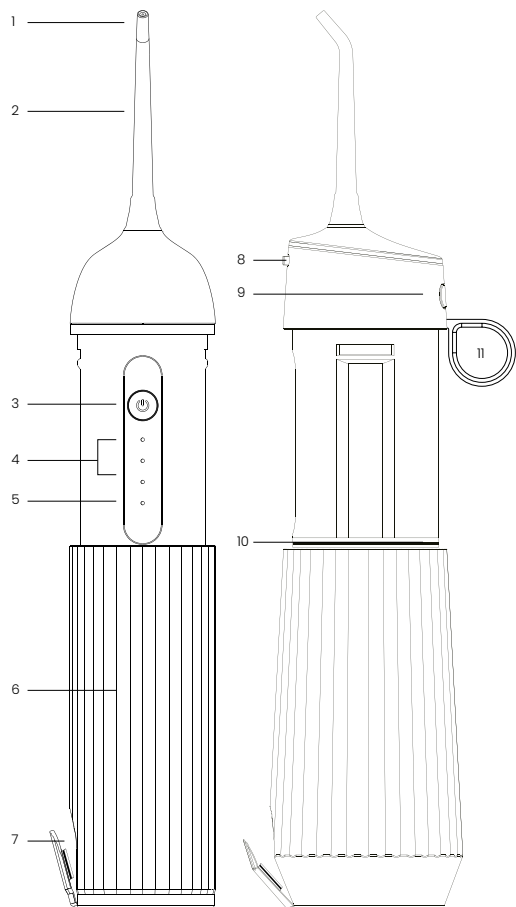PRODUCT GUIDE/SPECS

- ① Water Outlet
- ② CB-150H/150B nozzle
- ③ Power button
- ④ Mode indicator light
- ⑤ Battery indicator
- ⑥ Water storage tank
- ⑦ Storage tank cover & closed ring
- ⑧ Nozzle release button
- ⑨ Type-C charging port
- ⑩ Airtight ring of fuselage
- ⑪ Soft rubber handle

HOW TO CHARGE THE WATER FLOSSER

1. Connecting the Charger:

- Connect the Type-C end of the cable to the water flosser.
- Plug the USB end into a 5V power adapter (mobile phone adapter or computer).

2. Charging Tips:

- Charge for over 5 hours before first use or if not used for over 3 months.
- A flashing battery indicator means low battery.
- Ensure the charging cable is dry before use.
- Do not immerse the charging cable in liquid.

3. After Charging:

- The battery indicator will light up when fully charged.
- In soft mode, the device can be used for about 90 minutes.

HOW TO USE THE WATER FLOSSER

1. Preparing the Device:

- Insert the nozzle into the top of the product until you hear a click.
- Open the water tank lid and fill it with water.
- Close the water tank lid.

2. Using the Device:

- Place the nozzle above the gum line between teeth.
- Slightly close your lips to avoid splashing.
- Press the power/mode button to start the device (Soft, Standard or Pulse mode).
- Move the nozzle along the gum line and between teeth.

TROUBLESHOOTING

If the device does not work:

- Charge the device for 3 hours.
- Contact customer service if it still doesn't work.

If water pressure is low:

- Replace the nozzle.
- Ensure the water tank is filled.
- Use the device in a vertical position.

IMPORTANT SAFETY INSTRUCTIONS

When not to use the Product:

- If the power cord or plug is damaged.
- If the product does not work properly or is broken.
- If the product has been dropped into liquid.

User Restrictions:

- Not recommended for individuals (including infants and young children) with physical or sensory limitations, unless supervised.
- Keep out of reach of children and do not allow them to play with it.

Precautions to Avoid Injury:

- Do not aim the water jet under the tongue, in the ears, nose, or eyes.
- Consult a dentist if your gums bleed heavily after use or continue to bleed after 4 weeks.
- Consult a dentist if you have had recent oral or gum surgery or other dental treatments.

Medical Conditions:

Consult a doctor before use if you have medical conditions or use a pacemaker or other implanted device.

General Safety:

- Use only accessories recommended by Lucent.
- Do not disassemble the product.
- Do not pour hot water or unapproved solutions into the water tank.
- Do not insert foreign objects into openings.
- Do not clean with a dishwasher or microwave.

This product is for personal use only.

CLEANING AND MAINTENANCE

1. Cleaning the Nozzle:

- Press and hold the nozzle release button to remove the nozzle.
- Rinse the nozzle after each use.

2. Cleaning the Device:

- Remove the main unit from the water tank shell.
- Use a soft cleaning tool for the water tank and the device surface.
- If not using the product for a long period, clean and store it in a dry place.

3. Battery Disposal:

- Follow local regulations for battery disposal.
- Do not attempt to replace the battery yourself.

WARRANTY INFORMATION

All Be Lucent products are designed and manufactured to meet high expectations and to provide the best customer experience. However, sometimes, products may have malfunctions due to manufacturing defects. That is why we offer a full one year warranty from the date of original purchase, which covers defects in materials and workmanship.

The replacement warranty is valid only if the original receipt issued to the customer, specifying the date of purchase, is presented with the product to be replaced.

If you believe this product is defective, carefully repack the unit and return the product to your authorized seller within the warranty period, with the original receipt. Your authorized seller will replace the product, provided a defect in materials or workmanship is identified.

This limited replacement warranty does not apply in cases of abuse or misuse of the product, use contrary to the Lucent instructions, ordinary wear and tear, incorrect connection, force majeure or unauthorized repair.

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