©2026 STYLIDEAS Ltd. All rights reserved. STYLIDEAS owns trademarks, registered trademarks, patents, patents pending, design registrations and protection of the copyright of its products. www.stylideas.com/ip. Distributed by STYLIDEAS Ltd, UK: Amshold House, IG10 2RW, UK. +44(0)20 3225 5559 teamtom@stylideas.com. EU: OBELIS S.A, Bd. Général Wahis 53, B-1030 Brussels, Belgium

MC01B Made in China.

STYLPRO™ MINI MICROCURRENT Facial Toning Device

Appareil de Tonification Faciale par Microcourants
Microcurrent gezichtstonend apparaat.
Dispositivo per tonificazione viso a microcorrente
Urządzenie do Tonizacji Twarzy Mikroprądem
Dispositivo de Tonificación Facial con Microcorriente
Mikrostrøm Ansiktsoppstrammingsenhet
Dispositivo de Tonificação Facial por Microcorrente
Mikrovirta Kasvojen Kiinteytyslaite
Mikrostrøm Ansigtsoptæmmingsenhet
Mikrostrøm Ansiktslyftningsapparat

(EN) USE & CARE GUIDE

FOR YOUR SAFETY AND CONTINUED ENJOYMENT OF THIS PRODUCT, ALWAYS READ INSTRUCTION LEAFLET CAREFULLY BEFORE USING.

The product is intended for private use only and is not intended for commercial use or for use in other applications. This is not a medical product or a medical device.

SPECIFICATIONS


Name: STYLPRO™ MINI MICROCURRENT
Model: MC01B
Device Size: 57mm x 52mm x 24mm
Materials: ABS Plastic
Nominal Voltage: 3.7V
Rated capacity of the battery: 300mAh
Energy Capacity of the battery: 1.11Wh
Charge Voltage: 4.2V
Input power supply: USB socket, 5V, 1A
Charging Time: 2 hours

INCLUDES:

- 1 x STYLPRO Mini Microcurrent Device
- 1 USB-C Cable

PARTS IDENTIFICATION:

1. Microcurrent Device 2. On/Off Button
3. Intensity Level Indicators 4. USB-C Port
5. USB-C Charging Cable.

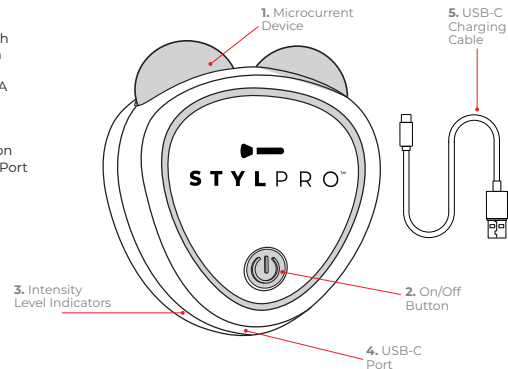3.7 V , 300mAh, 1.11Wh

SAFETY GUIDELINES:

- 1. We do not advise using this device if you are epileptic, if you are pregnant or have had skin cancer in the last 12 months, have undergone facial plastic surgery, have a metal implant, or if you are currently on your menstrual cycle.
- 2. For adult use only. Not intended for use by children. This is not a toy.
- 3. Wipe gently with a damp cloth.
- 4. Avoid storing the STYLPRO™ MICROCURRENT in hot, humid, or damp areas.
- 5. If the STYLPRO™ MICROCURRENT malfunctions or is damaged in any way, stop using it immediately. Contact STYLPRO™ CUSTOMER SERVICE (See section LIMITED TWO-YEAR WARRANTY for contact details).
- 6. Do not use on sensitive, irritated or broken skin. Do not use on swollen or sunburnt skin or if you have any other skin conditions.

CHARGING:

INSERT THE USB-C CABLE (PROVIDED) INTO THE PORT AT THE BOTTOM OF THE DEVICE. ATTACH TO A USB FRIENDLY POWER SOURCE. THE RED LIGHTS ON THE SIDE OF THE DEVICE WILL GLOW WHEN CHARGING.

HOW TO USE:

REMOVAL OF ANY MAKE-UP IS ADVISABLE BEFORE USE

- 1. ENSURE YOUR FACE AND NECK IS MAKEUP FREE AND CLEANSED.
- 2. APPLY YOUR FAVOURITE SERUM ON YOUR FACE AND NECK. THE MICROCURRENTS ARE LESS EFFECTIVE ON DRY SKIN. DO NOT USE THE DEVICE WITHOUT SERUM.
- 3. LONG PRESS THE ON/OFF BUTTON TO TURN ON THE DEVICE. YOUR DEVICE WILL BEEP, YOU WILL SEE 1 WHITE LED LIGHT TURN ON AT THE SIDE OF THE DEVICE AND A RED GLOW FROM THE CENTRE OF THE DEVICE.
- 4. TO ADJUST THE MICROCURRENT INTENSITY, SHORT PRESS THE ON/OFF BUTTON. THERE ARE 3 MICROCURRENT SETTINGS. ALWAYS START ON INTENSITY SETTING 1.
- 5. THE WHITE LED LIGHTS ON THE SIDE OF THE DEVICE WILL INDICATE WHICH MICROCURRENT SETTING YOU ARE CURRENTLY ON.
- 6. PRESS BOTH METAL SPHERES ON TO YOUR SKIN AND SLOWLY GLIDE THE DEVICE UPWARDS. GLIDE ACROSS YOUR CHEEKBONES, FOREHEAD, JAWLINE AND NECK. THE FIRMER YOU PRESS THE DEVICE, THE MORE POWERFUL THE MICROCURRENTS WILL FEEL.
- 7. GENTLY GLIDE OVER DESIRED AREAS FOR A WHOLE 10 MINUTE TREATMENT.
- 8. LONG PRESS THE ON/OFF BUTTON TO TURN OFF THE DEVICE. WIPE CLEAN WITH A DAMP CLOTH.

CLEANING:

WIPE WITH A DAMP CLOTH AND ALLOW TO DRY AFTER EVERY USE.

SAFETY INFORMATION :

Important!

Our products are low-risk, general wellness / fitness products that do not require medical clearance. They are designed as "General Wellness" home use products only.
After using the device, you may apply a moisturizer or other skincare products as part of your routine.

SAFETY CONSIDERATIONS:

Use of this device in any way other than those described in this manual voids the warranty and could result in serious injury.
The STYLPRO™ MINI MICROCURRENT is an electrical device. To avoid electric shock, adhere to the following safety considerations:

- Do not get your device wet or submerge your device in water.
- Do not wash electrical parts with water or liquids.

- Do not directly touch the device or switches with wet hands.
- Do not pull, carry, or lift the device by its charging cord. If the charging cord is damaged, do not use the device.
- Do not use the device if it has been dropped, left outside, or has encountered water or liquids.
- Do not use on or near heated surfaces.
- Do not modify, or attempt to modify or repair, this equipment. There are no user serviceable parts.
- Do not use or store your device in hot temperatures that are or can exceed 54°C (130°F).

IMPORTANT SAFETY INSTRUCTIONS: WHEN USING ELECTRICAL APPLIANCES, BASIC SAFETY PRECAUTIONS SHOULD ALWAYS BE TAKEN, INCLUDING THE FOLLOWING: READ ALL INSTRUCTIONS BEFORE USING APPLIANCE. SAVE THESE INSTRUCTIONS. KEEP AWAY FROM WATER. THIS PRODUCT IS NOT A TOY. USE INDOORS ONLY. FOR DOMESTIC HOUSEHOLD USE ONLY. DO NOT ATTEMPT TO OPEN THE BACK OF THE DEVICE. IF REQUIRED, CONSULT A PROFESSIONAL. KEEP AWAY FROM SOURCES OF DIRECT HEAT AND SUNLIGHT. KEEP AWAY FROM CHILDREN.

DANGER: TO REDUCE THE RISK OF ELECTRIC SHOCK: ALWAYS TURN OFF BEFORE CLEANING OR STORING.

WARNING:

TO REDUCE THE RISK OF BURNS, ELECTROCUTION, FIRE, OR INJURY TO PERSONS:

1. READ ALL INSTRUCTIONS BEFORE USING THE APPLIANCE.
2. THIS PRODUCT SHOULD NOT BE USED BY; CHILDREN, PREGNANT WOMEN, PATIENTS WITH EPILEPSY, TUMORS AND CANCER PATIENTS, PATIENTS WITH CARDIAC PACEMAKERS OR SIMILAR.
3. REMOVE GLASSES, FALSE EYELASHES AND CONTACT LENSES BEFORE USE.
4. IF YOU FEEL ABNORMALITIES IN YOUR SKIN DURING USE OR AFTER USE, PLEASE CONSULT YOUR MEDICAL SPECIALIST.
5. TO PROTECT AGAINST ELECTRICAL HAZARDS, DO NOT USE WHILE BATHING OR IN THE SHOWER.
6. CLOSE SUPERVISION IS NECESSARY WHEN ANY APPLIANCE IS USED BY OR NEAR CHILDREN OR PERSONS MADE WEAK OR DISABLED BY ILLNESS OR INJURY.
7. DO NOT USE ON SENSITIVE OR BROKEN SKIN. DO NOT USE ON SWOLLEN OR SUNBURNT SKIN OR IF YOU HAVE ANY OTHER SKIN CONDITIONS.
8. DO NOT OPERATE THIS APPLIANCE IF THE APPLIANCE MALFUNCTIONS OR HAS BEEN DROPPED OR DAMAGED IN ANY MANNER. RETURN APPLIANCE TO NEAREST SERVICE FACILITY FOR EXAMINATION, ELECTRICAL OR MECHANICAL ADJUSTMENT, OR REPAIR. NO REPAIRS SHOULD BE ATTEMPTED BY THE CONSUMER.
9. USE APPLIANCE ONLY FOR INTENDED USE AS DESCRIBED IN THIS INSTRUCTION BOOK. THE USE OF ATTACHMENTS NOT RECOMMENDED BY STYLPRO™ MAY CAUSE HAZARDS.
10. DO NOT USE OUTDOORS OR OPERATE WHERE AEROSOL (SPRAY) PRODUCTS ARE BEING USED, OR WHERE OXYGEN IS BEING ADMINISTERED.
11. DO NOT PLACE ANY PART OF THE PRODUCT NEAR HEAT SOURCES SUCH AS RADIATORS, COOKERS, ETC DO NOT LAY ON A SOFT SURFACE SUCH AS A BED AS DEVICE MAY OVERHEAT.
12. KEEP THE DEVICE AWAY FROM CORROSIVE LIQUIDS.
13. ENSURE DEVICE IS 'OFF' BEFORE STORAGE.
14. DO NOT USE APPLIANCE IN THE BATHROOM OR NEAR WATER.
15. DO NOT PLACE OR STORE APPLIANCE WHERE IT CAN FALL OR BE PULLED INTO A TUB OR SINK.
16. DO NOT PLACE IN OR DROP INTO WATER OR OTHER LIQUIDS.
17. IF THE APPLIANCE FALLS INTO WATER, DO NOT REACH INTO THE WATER.
18. CHOKING HAZARD DUE TO SMALL PARTS. NOT SUITABLE FOR CHILDREN UNDER 3 YEARS.
19. DO NOT SHARE THIS PRODUCT WITH ANYONE ELSE.
20. PLEASE DO NOT USE DURING CHARGING.
21. INSPECT THE DEVICE FOR ANY DAMAGE BEFORE EACH USE.

22. CONSULT WITH A DERMATOLOGIST OR HEALTHCARE PROVIDER IF YOU HAVE SENSITIVE SKIN, A HISTORY OF SKIN CONDITIONS, OR ARE UNCERTAIN ABOUT THE SUITABILITY OF THIS PRODUCT FOR YOUR SKIN TYPE.

BATTERIES: THE PRODUCT CONTAINS LITHIUM-ION BATTERY. RISK OF FIRE AND BURNS, DO NOT OPEN, CRUSH, HEAT ABOVE 60°C/140 °F OR INCINERATE. DO NOT SHORT CIRCUIT. IF THE LITHIUM-ION BATTERY BULGES SEVERELY, DISCONTINUE USE. KEEP SMALL CELLS WHICH ARE CONSIDERED SWALLOWABLE OUT OF REACH OF CHILDREN. SWALLOWING MAY LEAD TO BURNS, PERFORATION OF SOFT TISSUE, AND DEATH. SEVERE BURNS CAN OCCUR WITHIN 2 HOURS OF INGESTION. IN CASE OF INGESTION OF A CELL, SEEK MEDICAL ATTENTION PROMPTLY. NEVER RECHARGE NON-RECHARGEABLE BATTERIES. THIS PRODUCT HAS A BUILT-IN RECHARGEABLE BATTERY WHICH CANNOT BE REPLACED BY THE USER. THE REMOVAL OR REPLACEMENT OF THE RECHARGEABLE BATTERY MAY ONLY BE CARRIED OUT BY THE MANUFACTURER OR HIS CUSTOMER SERVICE OR BY A SIMILARLY QUALIFIED PERSON IN ORDER TO AVOID HAZARDS. WHEN DISPOSING OF THE PRODUCT, IT SHOULD BE NOTED THAT THIS PRODUCT CONTAINS A RECHARGEABLE BATTERY.

STORAGE

ENSURE DEVICE IS 'OFF' BEFORE STORAGE. ENSURE THE DEVICE IS UNPLUGGED FROM ANY POWER SUPPLY. PLACE IN COOL, DRY AREA OUT OF REACH OF CHILDREN, WHERE IT WILL NOT BE DAMAGED.

MAINTENANCE

YOUR STYLPRO™ MINI MICROCURRENT IS VIRTUALLY MAINTENANCE FREE. TO CLEAN TURN THE APPLIANCE OFF AND WIPE EXTERIOR WITH A DAMP CLOTH. IF ANY ABNORMAL CONDITION OCCURS, TURN OFF & CONTACT STYLPRO™ CUSTOMER SERVICES ON +44 (0) 20 3225 5559.

LIMITED TWO-YEAR WARRANTY

OUR GOODS COME WITH GUARANTEES THAT CANNOT BE EXCLUDED UNDER THE EUROPEAN CONSUMER LAW. EU DIRECTIVE 1999/44/EC, THIS DIRECTIVE STATES 'A TWO-YEAR GUARANTEE APPLIES FOR THE SALE OF ALL CONSUMER GOODS EVERYWHERE IN THE EU. YOU ARE ENTITLED TO A REPLACEMENT OR REFUND FOR A MAJOR FAILURE YOU ARE ALSO ENTITLED TO HAVE THE GOODS REPAIRED OR REPLACED IF THE GOODS FAIL TO BE OF ACCEPTABLE QUALITY AND THE FAILURE DOES NOT AMOUNT TO A MAJOR FAILURE. STYLPRO™ WILL REPAIR OR REPLACE (AT OUR OPTION) YOUR APPLIANCE FREE OF CHARGE FOR TWENTY-FOUR MONTHS FROM THE DATE OF PURCHASE, IF THE APPLIANCE IS DEFECTIVE IN WORKMANSHIP OR MATERIALS. THIS WARRANTY IS AN ADDITIONAL BENEFIT AND DOES NOT AFFECT YOUR LEGAL RIGHTS. TO OBTAIN SERVICE UNDER THIS WARRANTY, PLEASE CONTACT STYLPRO™ CUSTOMER SERVICE AT TEAMTOM@STYLIDEAS.COM +44 (0) 20 3225 5559. TO MAKE A WARRANTY CLAIM YOU MUST RETAIN YOUR PROOF OF PURCHASE.

IMPORTANT: THIS WARRANTY DOES NOT COVER DAMAGE RESULTING FROM ACCIDENTS, MISUSE, ALTERATIONS, ABUSE, LACK OF REASONABLE CARE OR USE OTHER THAN DESCRIBED IN THESE INSTRUCTIONS, OR PERFORMANCE OF SERVICE OR REPAIRS BY AN UNAUTHORISED SERVICE CENTRE. THIS WARRANTY IS VALID ONLY ON MERCHANDISE SOLD BY AUTHORISED RETAILERS. NOTE, THE ITEM MUST BE THOROUGHLY CLEANED PRIOR TO RETURN.

CORRECT DISPOSAL OF THIS PRODUCT: WHEN THE PRODUCT HAS REACHED THE END OF ITS USAGE, PLEASE DO NOT DISPOSE WITH REGULAR HOUSEHOLD WASTE. RECYCLING IS THE SAFEST AND MOST RESPONSIBLE WAY TO DISPOSE OF ELECTRICAL ITEMS THAT COME TO THE END OF THEIR LIFECYCLE. PLEASE DISPOSE OF IT ACCORDING TO YOUR LOCAL REGULATIONS. THE WASTE ELECTRICAL AND ELECTRONIC EQUIPMENT (WEEE) DIRECTIVE IN EUROPEAN COUNTRIES STIPULATES CORRECT RECYCLING PRACTICES. DO NOT DISPOSE IN REGULAR HOUSEHOLD WASTE.

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